VAC EN

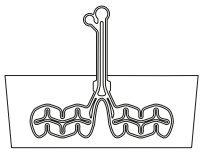
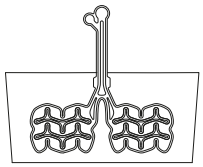
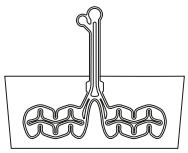
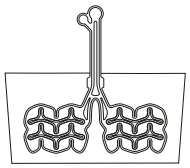
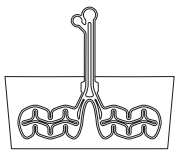
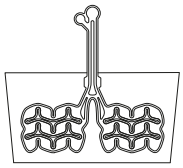
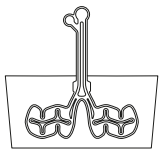
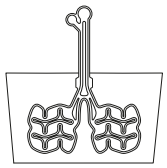
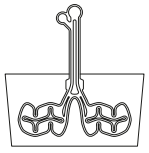
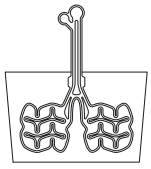
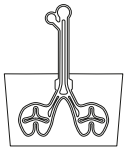
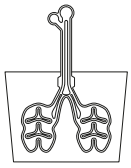
OVERLAY TF

SIZE TABLE

Smart Size Guide

Click here

This table provides an overview of all Overlay TF VAC sizes and their product dimensions. Please refer to the following pages for the patient measurement protocols.

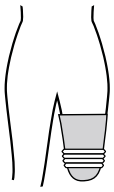
OVERLAY TF CIRCUMFERENCE* *Product circumference at 4 cm from bottom of sleeve	SH 20 cm	LG 25 cm
52 52 cm		
48 48 cm		
44 44.5 cm		
40 40 cm		
38 36 cm		
32 32 cm		

OVERLAY TF LENGTH*

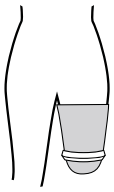
*Length from bottom of sleeve to top

INTENDED FOR :

Single seal



Multi seal



ETHNOCARE™

Please see next page for **PIN** table

PIN EN

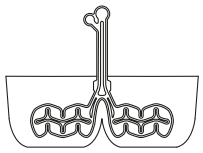
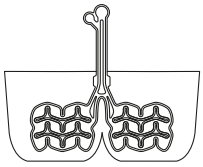
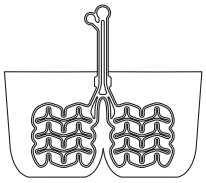
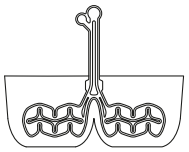
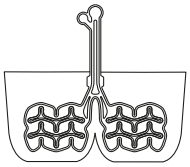
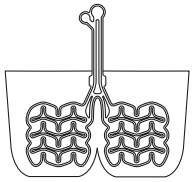
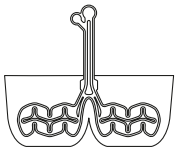
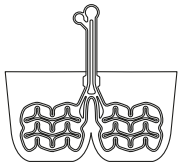
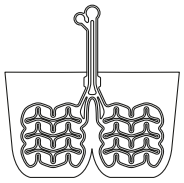
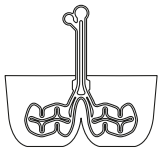
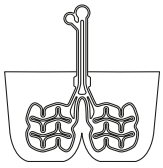
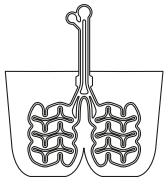
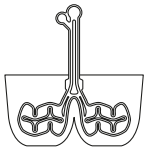
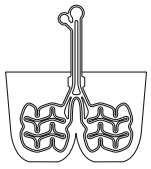
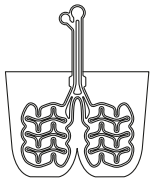
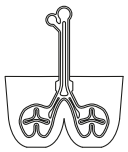
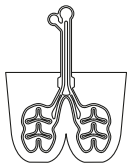
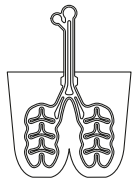
OVERLAY TF

SIZE TABLE

Smart Size Guide

Click here

This table provides an overview of all Overlay TF PIN sizes and their product dimensions. Please refer to the following pages for the patient measurement protocols.

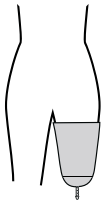
OVERLAY CIRCUMFERENCE* *Product circumference at 4 cm from bottom of sleeve	SH 20 cm	LG 25 cm	XL 30 cm
52 54 cm			
48 50 cm			
44 47 cm			
40 42 cm			
38 38 cm			
32 28 cm			

OVERLAY TF LENGTH*

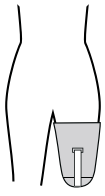
*Length from bottom of sleeve to top

INTENDED FOR :

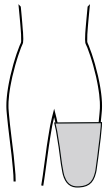
Pin



Lanyard



Cushion



ETHNOCARE™

Please see next page for **PIN** sizing

PIN **EN**

Compatible with pin, lanyard and cushion liners.

OVERLAY TF

SIZING CHART

To determine the correct size, measure the limb's circumference and length at a specific height.

Smart Size Guide

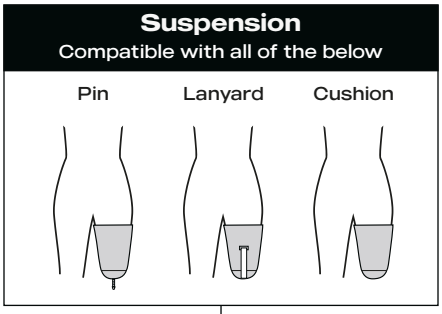
Click here

STEP 1

Choose liner type.

Suspension

Compatible with all of the below



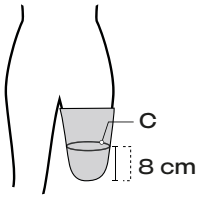
PIN

STEP 2

Measure circumference.

Circumference

Measure at **8** cm from the distal end



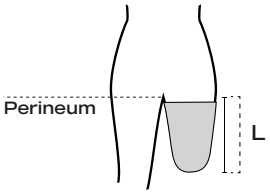
Measurement	32-36 cm	36-40 cm	40-44.5 cm	44.5-48 cm	48-52 cm	52-62 cm
SKU	32	38	40	44	48	52

STEP 3

Measure length.

Length

Measure from the **distal end** of the liner to the **perineum**



For check sockets you can also measure to the medial trim line.

Measurement	20-25 cm	25-30 cm	30+ cm
Bladder rows	2	3	4
SKU	SH	LG	XL

Suspension	Circumference	Length	SKU
PIN	32-36 cm	20-25 cm	OVTF-32-SH-PIN
	36-40 cm		OVTF-38-SH-PIN
	40-44.5 cm		OVTF-40-SH-PIN
	44.5-48 cm		OVTF-44-SH-PIN
	48-52 cm		OVTF-48-SH-PIN
	52-62 cm		OVTF-52-SH-PIN
	32-36 cm	25-30 cm	OVTF-32-LG-PIN
	36-40 cm		OVTF-38-LG-PIN
	40-44.5 cm		OVTF-40-LG-PIN
	44.5-48 cm		OVTF-44-LG-PIN
	48-52 cm		OVTF-48-LG-PIN
	52-62 cm		OVTF-52-LG-PIN

Suspension	Circumference	Length	SKU
PIN	32-36 cm	30+ cm	OVTF-32-XL-PIN
	36-40 cm		OVTF-38-XL-PIN
	40-44.5 cm		OVTF-40-XL-PIN
	44.5-48 cm		OVTF-44-XL-PIN
	48-52 cm		OVTF-48-XL-PIN
	52-62 cm		OVTF-52-XL-PIN

Please see next page for **VAC** sizing

VAC EN

Compatible with vacuum and suction suspension systems using a sealing liner.

OVERLAY TF

SIZING CHART

To determine the correct size, measure the limb's circumference and length at a specific location depending on sealing liner type.

Smart Size Guide

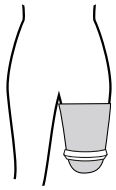
Click here

STEP 1

Choose liner type.

Single seal

1 seal, sealing liners



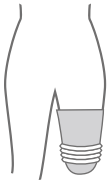
Overlay will sit under the sealing ring

VAC

Multi seal

4 or 5 seal, sealing liners

Follow instructions for single seal but measure **circumference (C)** just above top seal and **length (L)** from that location to perineum



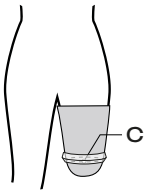
Overlay will sit above the top sealing ring

STEP 2

Measure circumference.

Circumference

Measure **under** the reflected sealing ring



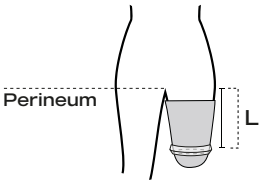
Measurement	32-38 cm	38-42 cm	42-47 cm	47-50 cm	50-54 cm	54-62 cm
SKU	32	38	40	44	48	52

STEP 3

Measure length.

Length

Measure **under** the reflected sealing ring to the perineum



Measurement	20-25 cm	25+ cm
Bladder rows	2	3
SKU	SH	LG

Suspension	Circumference	Length	SKU
VAC	32-38 cm	20-25 cm	OVTF-32-SH-VAC
	38-42 cm		OVTF-38-SH-VAC
	42-47 cm		OVTF-40-SH-VAC
	47-50 cm		OVTF-44-SH-VAC
	50-54 cm		OVTF-48-SH-VAC
	54-62 cm		OVTF-52-SH-VAC

Suspension	Circumference	Length	SKU
VAC	32-38 cm	25+ cm	OVTF-32-LG-VAC
	38-42 cm		OVTF-38-LG-VAC
	42-47 cm		OVTF-40-LG-VAC
	47-50 cm		OVTF-44-LG-VAC
	50-54 cm		OVTF-48-LG-VAC
	54-62 cm		OVTF-52-LG-VAC